



TX® Traction System

4759 – TX Traction Unit

Standard Accessories

48082 – Patient Interrupt Switch

27465 – Patient Data Cards (5x)



New innovation for a classic therapy

In 1979 Chattanooga Group introduced the TX line of clinical traction devices to the physical medicine market. Today, Chattanooga Group is the world leader in clinical traction devices, tables and accessories. Using the most advanced technology, the new TX traction unit is years ahead of any traction unit currently available on the market today. With innovative features and an easy to use digital touch screen interface that makes set up simple, the redesigned TX traction unit is the first true innovation to this classic therapy in over 25 years.


TX Innovative Features

- Easy to operate
- Digital Touch Screen Interface
- 270° Pivoting Interface
- Patient Pain Mapping
- Patient Data Card records pre- and post-treatment pain profiles, treatment session information and patient information
- 10 user-defined protocols
- Real-time digital treatment tracking
- User-defined Hold, Rest and Treatment times
- Intermittent, Static and Cyclic traction programs
- Progressive and Regressive settings
- Traction tension parameters 0 – 200 lbs.
 - » Patient Interrupt Switch
 - » Cervical maximum poundage warning
 - » Audible signal at end of treatment or when Patient Interrupt Switch is activated

USER PROTOCOLS



CLINICAL LIBRARY

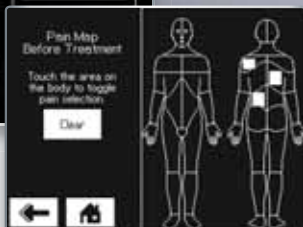


PAIN MAP

Technical Specifications

Mains Power: 100 – 240 V, 50/60 Hz
 Power Consumption: 75 W
 Current Consumption: 3.2 Amps Max
 Weight: 30 lbs (14 kg)
 Dimensions: 17.5" L x 9.5" W x 17.5" H
 (45 cm x 24 cm x 45 cm)
 Electrical Safety Class: Class 1, Type B
 Safety Tests: IEC/UL/EN: 60661-1, 60601-1-2;
 CAN C22.2 No. 601.1-M90 w/AZ;
 cETLus Listed



TRACTION TENSION PARAMETERS	MIN - MAX	INCREMENTS
TRACTION PERIOD	1 - 99 MINUTES	1 MINUTE
HOLD PERIOD	0 - 99 SECONDS	1 SECOND
REST PERIOD	0 - 99 SECONDS	1 SECOND
TRACTION TENSION	0 - 200 LB (0 kg/0 N) - (90kg/890 N)	1 LB (1 kg/5 N)
PROGRESSIVE STEPS	1 STEP - 9 STEPS	1 STEP
REGRESSIVE STEPS	1 STEP - 9 STEPS	1 STEP

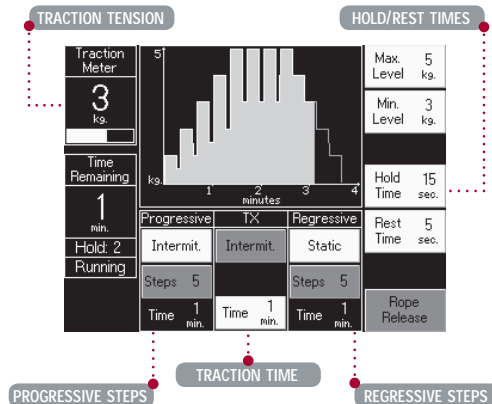

Patient Data Management System

2790 – Patient Data Management System (PDMS)
 Includes: software, card reader and 5 patient data cards

27465 – Patient Data Cards (25/pack)

Mobile Traction Stand

3600 – Mobile Traction Stand

Patient Data Management System (PDMS)



PDMS is a unique optional software program that allows the clinician to organize, store, display and print patient traction therapy information from the TX device to the PC. This system which includes software and a card reader for your PC can be used with the patient data cards (5 included with PDMS) to document treatments.

Mobile Traction Stand

- Fully adjustable
- Allows the TX Traction device to be used with most tables
- Traction unit sold separately

