



ENTHERMICS Medical Systems

EC1540bl

Combination Blanket/Fluid Warming Cabinet



Specifications:

Dual-cabinet design with 20 gauge stainless steel exterior casing and doors. Window doors and interior lighting allow observation of inventory with the doors closed and convenient illumination when the doors are open. The cabinet is warmed using the patented Halo Heat® low-heat-density electrothermal cable array. Fluids compartment is equipped with three (3) white, epoxy-coated wire baskets, each with a 14 liter capacity to accommodate most commonly used surgical fluids packaged in bags or bottles. Blanket compartment includes a white, epoxy-coated blanket support assembly. The cabinet is furnished with a full-perimeter rubber bumper assembly and one (1) set of 5" (127mm) heavy duty casters — two with locking brake.

Special Information:

- Right-hand door swing is standard.
Left-hand door is available as a special order.

Control:

Each compartment is controlled by one (1) electronic control which consists of a 4 digit L.E.D. display, ON/OFF key, heat indicator L.E.D., INCREASE and DECREASE keys, built-in lock out feature and a series of prompt sequence indicators. Controls can easily be set to operate in Celsius or Fahrenheit. In the event of a power failure, the cabinet will remember its programming and begin to operate as before. A warming shut-off system, separate from the electronic control, prevents overheating. The fluids compartment has individual temperature controls for irrigation fluids (IRR) and injection fluids (INJ) to allow for precise temperature settings. IRR temperature may be adjusted from 90° to 150°F (32° to 66°C) while the INJ temperature range is 90° to 110°F (32° to 43°C). An alarm will sound if temperatures exceed 10°F (6°C) above the set-point temperature. The blanket compartment thermostat has an adjustable temperature range of 90° to 200°F (32° to 94°C). Auto mode allows the user to program the control to automatically turn on and turn off once during a 24 hour period at a selected time.

Capacity and Weights


Interior Capacity: (77 cubic feet per compartment)	
Fluid Compartment	42 liters
Blanket Compartment	26-5/8" x 18-7/8" x 26-7/8" (676mm x 479mm x 683mm)
Net Weight:	412 lb (187 kg)
Ship Weight:	502 lb (228 kg)

Electrical

120 V.A.C. 50/60 Hz, 1ph
1920 Watts, 16.0 Amps



230 V.A.C. 50/60 Hz, 1ph
1700 Watts, 7.5 Amps

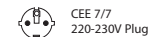

Options and Accessories

Combination Lock	5004663
Cylinder Lock for Door Handle	LK-22567
Legs Kit, 6" (152mm)	44093
Window Door Assembly, LH, ea.	5005863

