

Ultra Aire Select™

Bariatric Low Air Loss and Alternating Pressure Therapy



Model Number: P2500-39, 42, 48, 54

Ultra Aire Select™

Dynamic Low Air Loss and Alternating Pressure Therapy for Large or Bariatric Patients

Control unit is sleek, lightweight, and easily attaches to the footboard.



The Control Unit

- Capable of supporting patients weighing up to 1000 lbs.
- Provides effective low air loss to assist with moisture control.
- Selectable Alternating Pressure with a five minute cycle time.
- User-friendly controls allow the caregiver to adjust for individual patient comfort.
- Equipped with Max. Inflate for easy patient positioning.
- Patient settings are restored after power outages.
- CPR quick release for rapid deflation.
- Compact and lightweight.


The Therapeutic Support Surface

- All air cushions are constructed of high quality and very pliable coated nylon for low shear and even lower interface pressures.
- Plexus® Quilted Top Sheet is highly vapor permeable and oversized providing a low friction, low shear barrier against complications due to moisture.
- Unique modular design facilitates cleaning and interchangeability.
- Convoluted foam base provides additional safety for your patient.

FootStop/ZIP™ Cushion (Optional)

- Interchangeable Foot Stop/ZIP™ cushion permits tailoring the mattress length to fit the patient while reducing pressure in the critical heel area.
- ZIP cushion can be independently used anywhere along the mattress.

Optional Foot Stop/ZIP™ cushion uniquely maintains proper patient positioning while reducing interface pressure on the critical heel area.

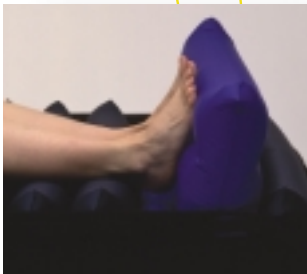

The therapeutic support surface is lightweight and made of high quality materials. It conveniently rolls up for easy carrying with attached straps. Separate carrying case also available.



Specifications: P2500-39, 42, 48, 54 Mattress Replacement System

Dimensions:	Control Unit: 5" x 10" x 12" / 12.7cm x 25.4cm x 30.5cm (9.5lbs. / 4.3kg) Mattress: 80" x 10" x (39", 42", 48" or 54" width) 203.2cm x 25.4cm x (99cm, 107cm, 122cm or 137cm width)
Electrical Rating:	120V~, 60Hz, 5A; / 230V~, 50Hz, 2.5A
Agency Approvals:	UL 544 C22.2 No. 125; TUV, EN 60601-1
AP Cycle Time:	5 minutes
Warranty:	* Control Unit: 2 year Ltd. Mattress: 1 year Ltd.
Weight Limit:	1000 lbs.
Medicare HCPCS Code:	E0277

*Available in 220 Volt model



GAYMAR®