



Lift Specifications

Dimensions	
Height	
Maximum at sling hookup	70"
Minimum at sling hookup	29.5"
Base height (clearance)	6.8"
Base width	
Open	41"
Closed	26.5"
Base length	48"
Casters	
Front	5"
Rear	5"
Weight capacity	600 lb.
Product weight	
Out of carton	112 lb.
Shipping	138 lb.

Sling Specifications



	R140	R141
Dimensions		
Width		
	45"	45"
Commode opening		
	n/a	7"
Length		
	55"	55"
Commode opening		
	n/a	13"
Product weight		
Out of carton	2.2 lb.	2.1 lb.
Shipping	2.7 lb.	2.6 lb.
Weight capacity	600 lb.	600 lb.


Invacare® Reliant 600 Heavy-Duty Lift

Model no. RPA600-1


Designed to meet the special needs of your patients weighing up to 600 pounds, the Invacare Reliant 600 heavy-duty lift is secure and sturdy, providing a safe alternative to manual lifting. Both durable and functional, the Reliant 600 is also a reassuring and comfortable patient lift. When answering the needs of every patient is your priority, this is your wisest choice.

FEATURES

- High lift range allows transfer to high surfaces
- Low friction casters significantly improve rollability and contributes to caregiver safety
- Designed to meet the needs of the growing bariatric population up to 600 lb.
- Manual lowering for instant response to power loss
- Safety tested with  mark for highest quality standards
- Essential device for safe patient handling in hospitals and LTC (long-term care) facilities
- Unique sling fabric designed to provide the utmost in comfort and strength
- Warranty includes three years on lift, one year on electric components
- Rear locking casters are standard
- Manufactured in an ISO9001 facility 

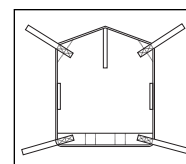

Model no. R140 Mesh Sling

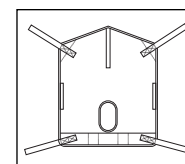

Model no. R141 Mesh Sling with Commode Opening

Complete instructions for lift usage and sling procedures are found in the Invacare Operating Manual. Only genuine Invacare slings and accessories should be used on Invacare lifts.




Invacare® Digital Scale

Model no. RSC600

- Compact precision scale mounted on swivel bar accurately measures weights up to 600 lb.
- Provides wide lift range and retrofits on all Reliant 450/600 lifts. Model no. 9805 requires Reliant swivel bar.



Charging Station

- One 24-volt battery on mast with quick release
- Battery can be replaced immediately with second battery in  approved, wall-mounted charger
- UL listed and CSA certified



Manual Lowering

Essential Safety Feature

- Manual emergency lowering feature gently and safely lowers patient in case of total power loss
- Manual lowering is easy and effortless with an accessible ring-pull

Invacare® Reliant slings come in a wide array of styles to meet a variety of client needs, and, the new, improved overall design makes them all easier than ever to position. Reliant slings feature a unique, fully padded fabric that forms a soft yet supportive surface. The soft-brushed lining helps prevent sliding for ultimate comfort and skin protection. And the new stretch-resistant feature means the slings will hold their shape, for superior support and ease of positioning. For clients and caregivers alike, Invacare Reliant slings are designed with comfort and safety in mind.



Invacare® Reliant Full-Body Slings

Model nos. R110, R111, R112, R113, R117

- Polyester/nylon fabric and/or mesh offers full head and neck support
- Four-point hook-up
- Fully padded
- Available in four sizes: M, L, XL, XXL



Invacare® Reliant Full-Body Slings with Commode Opening

Model nos. R114, R115, R116

- Mesh fabric offers full head and neck support
- Commode opening for toileting
- Four-point hook-up
- Available in four sizes: M, L, XL, XXL



Invacare® Reliant Divided-Leg Slings with Headrest

Model nos. R100P, R100, R101, R102

- New supportive back and contoured headrest are easily inserted or removed from wheelchair while patient is seated
- Excellent for all transfers, especially toileting and bed-to-chair
- Three support straps on each side for safety and security
- Fully padded
- Available in four sizes: P, M, L, XL



Invacare® Reliant Toileting Slings with Belt*

Model no. R121

- Easily removable firm back and contoured headrest
- Excellent for all transfers, especially toileting and bed-to-chair
- Three support straps on each side for safety and security



* Toileting sling allows superb access for hygiene and clothing removal, but patient must have full head and neck control for proper use of this sling.

Sling Specifications

	Divided-leg				Full-body			Full-body w/commode			Toileting
	P	M	L	XL	M	L	XL	M	L	XL	L
Dimensions											
Width	33"	37.5"	41"	44.5"	41.5"	45.5"	45.5"	41.5"	45.5"	45.5"	36"
Commode opening	n/a	n/a	n/a	n/a	n/a	n/a	n/a	8"	8"	9"	n/a
Length	59.8"	62.8"	67.8"	72.3"	54.7"	60.5"	65.3"	54.7"	60.5"	65.3"	37"
Commode opening	n/a	n/a	n/a	n/a	n/a	n/a	n/a	11"	11"	13"	n/a
Back	29.2"	35.7"	41.7"	43.2"	n/a	n/a	n/a	n/a	n/a	n/a	n/a
Product weight											
Out of carton	0.9 lb.	0.9 lb.	1.0 lb.	1.1 lb.	0.6 lb.	0.7 lb.	0.8 lb.	0.6 lb.	0.7 lb.	0.8 lb.	0.9 lb.
Shipping	1.5 lb.	1.5 lb.	1.6 lb.	1.7 lb.	1.2 lb.	1.3 lb.	1.4 lb.	1.2 lb.	1.3 lb.	1.4 lb.	1.5 lb.

Complete instructions for lift usage and sling procedures are found in the Invacare Operating Manual. Only genuine Invacare slings and accessories should be used on Invacare lifts.
 Note: Only Reliant sling models listed below are for use with lift Model nos. 9805, Reliant RPA450-I, RPL450-I, RHA450-I and RHL450-I.

Full-Body

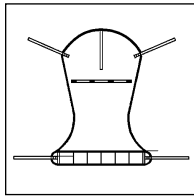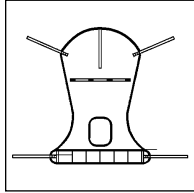

Model no.	Fabric / Style	Size
R110	mesh/full-body	M
R111	mesh/full-body	L
R112	solid/full-body	M
R113	solid/full-body	L
R114	mesh/full-body with commode opening	M
R115	mesh/full-body with commode opening	L
R116	mesh/full-body with commode opening	XL
R117	solid/full-body	XL

Note on heavy-duty and stand-up slings:
 R130 and R134 available for use on RPS350-I stand-up only. See page 4.
 R140 and R141 available for use on RPA600-I heavy-duty lift. See page 5.
 Standard slings for use with chains or straps. See back cover.


Sling Binding Colors

Refer to chart below.



Sling Sizes

Color / Weight Chart

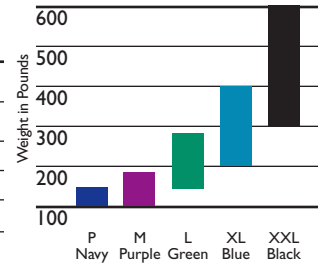


Note: Sling size and fit can vary significantly depending on weight and girth. These are general guidelines. Consult physician before sling selection.

Washing instructions are on each sling label. Washing temperature should not exceed 120°F. Air-dry or dry at low temperatures.


All slings are constructed of polyester/nylon fabrics which can be used in wet or dry environments.

Divided-Leg

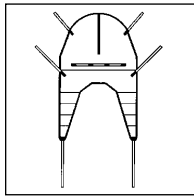

Model no.	Fabric / Style	Size
R100P	solid/divided-leg with headrest	P
R100	solid/divided-leg with headrest	M
R101	solid/divided-leg with headrest	L
R102	solid/divided-leg with headrest	XL

Toileting

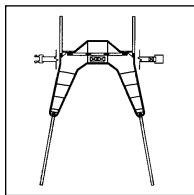

Model no.	Fabric / Style	Size
R121	solid/toileting	L

Note: Patients must have full head and neck control for proper use of this sling.

Swivel Bars

- Curved sling attachment safety hooks prevent straps from disengaging
- Accommodates all slings from webbing to chain-style
- Rotates 360° to increase maneuverability
- Design eliminates swing and sway
- 450-lb. or 600-lb. weight capacity
- Covered with protective padding (450 lb.) or protective coating (600 lb.)
- Accommodates digital scale, Model no. RSC600 (sold separately)
- Retrofit kit available for older lift models:
 - Part no. 1068992 (can be retrofitted to Model no. 9805)
 - Part no. 1075644 (can be retrofitted to Reliant 450 lifts only)



Standard Swivel Bar (Part no. 1068992)



ExtraWide Swivel Bar* (Part no. 1075644)

*Cannot retrofit to Model no. 9805.