


MODEL 1036-19 "L" SHAPED AUTOPSY TABLE




STANDARD DESIGN FEATURES:

- Autopsy area is 30 inches wide x 106 inches long
Dissecting area is 24 inches wide x 60 inches long
- All Type 304 stainless steel construction with a polished No. 4 Finish
- Hot and Cold Water Mixing Faucet with Gooseneck Spout and Wrist Blade Handles
- Deluxe Hydro-Aspirator with Reverse Flow, Cold Water Control Valve and 10 feet of hose
- Hand Spray Hose Assembly with Cold Water Control Valve, Nozzle and 10 feet of hose
- Table End Rinse Assembly with Cold Water Control Valve on autopsy section
- Heavy-Duty Commercial Waste Disposal (1/2 hp) with vacuum breaker and solenoid
- Large Single Compartment Sink
- Vacuum Breaker protection for all plumbing fixtures to prevent backflow of contaminated water
- Double Cone Down Draft Ventilation Assembly on autopsy section
- Full length perforated dissection grids for autopsy section
- Acoustical Coating to reduce metallic noise and prevent condensation
- Tabletop is sloped toward sink basin for rapid and positive drainage

UNITS ARE AVAILABLE AS LEFT HAND OR RIGHT HAND CONFIGURATIONS.

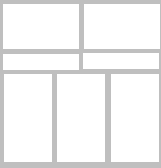

MODEL 1036-19 LAYOUT DRAWING and OPTIONAL FEATURES

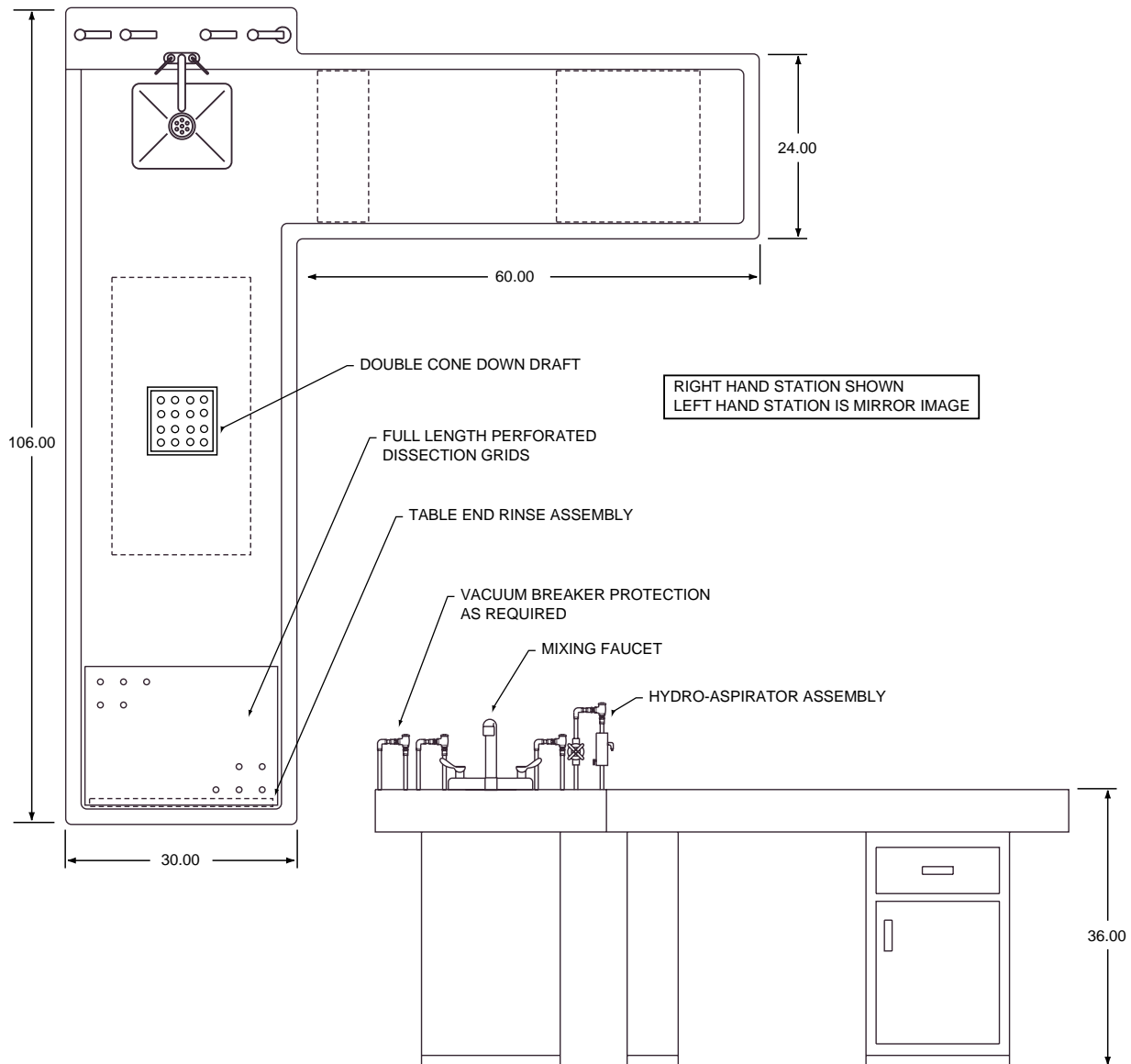

OPTIONAL EQUIPMENT AND ACCESSORIES:

- Heavy Duty Commercial Disposal (3/4, 1 or 2 hp)
- Post Mortem Scale with Scale Pole
- Gas / Air Cocks
- Emergency Eye Wash
- Centimeter or inch markings
- Over-the-body Dissecting Tray
- Headrest