



PressureGuard® APM²

Category: Alternating Pressure with Basic Lateral Rotation


Alternating Pressure or Lateral Rotation – at the Flip of a Switch!

With its two distinct therapy modes, the PressureGuard APM² gives you twice the treatment options. Twice the flexibility. Twice the value.

Change instantly and easily from alternating pressure to basic lateral rotation with a single, recessed toggle switch. Requires no disconnecting or exchanging of air lines.

New Deluxe control unit includes “static” button which parks the unit in powered flotation mode, temporarily suspending rotation/alternating pressure to facilitate feeding, transfers, dressing changes, etc. Low pressure indicator light gives users at-a-glance assurance of proper function.

The APM² delivers all the patented features that make our original PressureGuard APM stand out from the pack: Therapeutic Geo-Matt anti-shearing surface, two-part Safety Edge™, subtle shifts in weightbearing in ten-minute cycles without sharp pressure drop-offs.

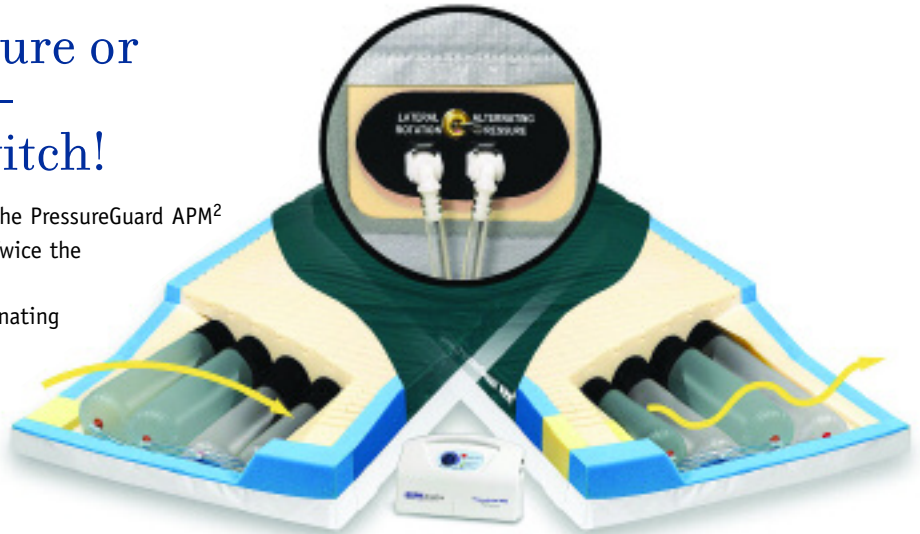

QUICK FACTS

PressureGuard® APM²/Safety Supreme

Primary Use:

- Aggressive Treatment of Pressure Ulcers

Patient Indications:

- Can't or won't reposition frequently
- Can be positioned off of ulcer in at least 2 positions
- Prevention for those at high risk for skin breakdown (e.g. Braden score 12 or less)

Key Advantages:

- Two therapies in one mattress
- Single switch operation
- PressureGuard® design: stability, safety, anti-shearing



Deluxe control unit

- “static” flotation mode
- low pressure indicator
- shielded quick-connects
- 10-minute cycle
- comfort control
- hospital-grade plug

APM² Models:

Item #	Description
5875LR-13	PressureGuard® APM ² , w/deluxe control unit, 75”L x 35”W x 7”H
5880LR-13	PressureGuard® APM ² , w/deluxe control unit, 80”L x 35”W x 7”H
5884LR-13	PressureGuard® APM ² , w/deluxe control unit, 84”L x 35”W x 7”H
5800	Deluxe control unit only
5500	Standard control unit only

Wider sizes available: See In-Home Styles – page 14

APM Models: Ships with standard control unit, which does not include lateral rotation feature or static mode.

Item #	Description
5175-13	PressureGuard® APM, w/standard control unit, 75”L x 35”W x 7”H
5180-13	PressureGuard® APM, w/standard control unit, 80”L x 35”W x 7”H
5184-13	PressureGuard® APM, w/standard control unit, 84”L x 35”W x 7”H

Bariatric APM — see next page

Medicare Approved: E0277 (Group II)
Weight limit: 350 lbs.
Warranty: 18 month, non pro-rated
Mattress Weight: 20 lbs.
UL listed: Yes

Fire barrier versions can be ordered to conform to NFPA 101 (Life Safety Code); Cal Tech Bulletins # 117, 129, and 603; 16 CFR Parts 1632 and 1633; and BFD IX-11. To order, replace “-13” suffix with “-29” suffix in model number.